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TEST REPORT

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NON-MEDICAL

How to Use

- 1. Open the mask, with the white side facing in and nose clip on the top.
- 2. Pull the ear straps behind your ears and adjust them to make you feel comfortable. 3. Use your hands to adjust the shape of the nose clip until it is pressed to fit the bridge of your nose.
- Adjust the size of the mask by pulling the top and bottom until you cover your nose and mouth.

Instruction: Cover on the nose and mouth

- Precautions:

 Before using the product, check whether the packaging is in good condition, confirm the production date of the packaging indication, and use within the validity period;

 This product is non-sterilized. If the package is damaged, please discontinue use; (If it is necessary to enter the sterile room or other environment, please enter after sterilization. Ethylene oxide is recommended for sterilization;

 This product is disposable and cannot be re-used. Please use as soon as possible after it is unsealed.

 The maximum wearing time should not exceed 4 hours;

 Please avoid direct contact with the inside of the mask during use;

 After use, the product should be disposed according to the provisions of the hospital or the environmental protection department;

 Please read the instruction carefully before use;

Caution: This product is not

Product Storage Conditions and Methods: This product should be stored in the room with relative humidity of less than 80%, good ventilation, no corrosive gas, fire, moisture, anti-rodent and anti-insect equipment, and the storage period should be within the product validity period;

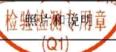
Valid Date: Storage under prescribed conditions, valid for two years, please use within the validity

See the packaging indication for production batch number and expiration date.
Valid Date: 2 Years

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Executive Standard: GB/T32610-2016

Specification: 175MM*95MM(+-5MM)

3 LAYERS ,50PCS/BOX(50PCS/BAG) LARGE SIZE

Product Structure and Main Purpose:

- The product of Earloop Mask is composed of the mask body, elastic band and nose clip.
- The mask body is composed of PP non-woven fabric, Melt-blown Fabric and PP non-woven fabric.

 The upper part of the mask is equipped with a nose clip made of flexible material.

Scope of Use:

it is suitable for the health care of the user in general environment or the open environment where there is no risk of body fluids and spillage this mask can not be use as a medical protective mask

MANUFACTURE: Dongguan Youai Medical Supplies Co., Ltd. FACTORY ADDRESS: 149-1 Qinglong Road, Huangjiang Town, Bongguan City, Guangdong Prevince, China

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